


# dreamTAP Cheat Sheet

## Key Features

- Single Point Adjustment
- Dental Grade Cr-Co Alloy
- Snap-Fit Trays
- 15mm protrusion range
- 15mm lateral excursion range
- Maximum lingual tongue space


## Protrusion Range

- 15 mm w/ three hooks

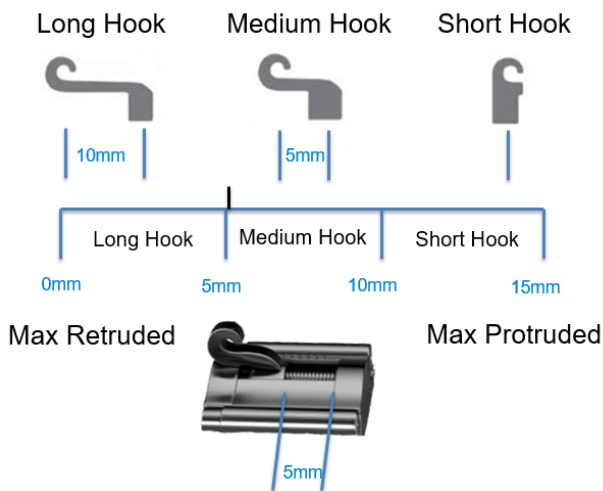


## Adjusting Protrusion

- Advance counter-clockwise
  - Facing patient
- Half turn = 0.25 mm
- More than 5mm, change hook

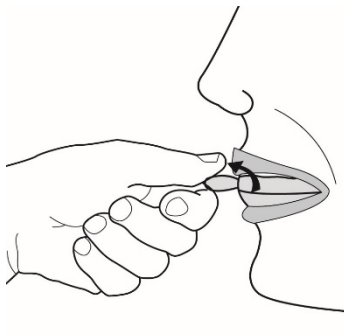

## Hardware Options

- Quick Release Bar
- Medicare Bar
- 4-15 mm Lateral Excursion

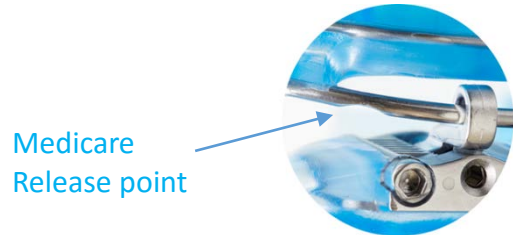

## Tray Options

- TL 1.8mm duo-durometer
- ThermAcryl 2 mm PMMA lined
- 5-15mm Vertical height
- Posterior Stops
- Heavy Bruxer Option
- Compliance Chip add-on



## Standard Design

- Quick release bar
- Set at 10mm forward & 5mm backward
- About 15mm lateral excursion
- About 5mm of vertical in the anterior
- TL Tray

<b>comparison chart</b>	<b>dreamTAP</b>	<b>Herbst</b>	<b>Dorsal Fin</b>	<b>Tray Style</b>	
Single Point Adjustment	<b>Yes</b>	<b>No</b>	<b>No</b>	<b>No</b>	Bi-lateral adjustment can cause TMD problems
Maximum Protrusion Range	<b>15 mm</b>	5-7 mm	5-8.5 mm	5-15 mm	Bigger range makes it easy to hit the protrusion target
Micro-Adjustments < .5 mm increments	<b>Yes</b>	<b>Yes</b>	<b>Yes</b>	<b>No</b>	Fine tuning makes it easy to hit the protrusion target
Lateral Excursion Range	<b>4-15 mm</b>	2-4 mm	0 mm	3-4 mm	Greater lateral excursion is more comfortable for patient
Vertical Range	<b>5-15 mm</b>				
Hardware Material	<b>Cr-Co</b>	SS	SS	Plastic	Dental grade Cr-Co is the best alloy for the oral environment, bar none
Tray Material	<b>Snap-Fit Polycarbonate</b>	Acrylic w/ball clasps	Acrylic w/ball clasps	Nylon/Plastic	Flexible trays create greater retention without teeth moving ball clasps
Prevents Mouth Breathing	<b>Yes</b>	<b>With rubber bands</b>	<b>No</b>	<b>No</b>	Eliminating mouth breathing improves efficacy
Seats 100% first time, every time	<b>Yes ThermAcryl</b>	<b>No</b>	<b>No</b>	<b>No</b>	ThermAcryl lining allows for a perfect chairside fit, every time
Compliance Chip Available	<b>Yes</b>	<b>Yes</b>	<b>Yes</b>	<b>Yes</b>	Future DOT requirement
Multiple Tray Options	<b>Yes</b>	<b>No</b>	<b>No</b>	<b>No</b>	A choice of trays to fit every patient
Add-on CPAP Mask for Combo Therapy	<b>Yes</b>	<b>No</b>	<b>No</b>	<b>No</b>	Add a mask to treat severe cases. Also improves MD relationships
Peer Reviewed Studies	<b>36</b>	2	2	2	
ISO 13485 Certified Manufacturing	<b>Yes</b>	<b>No</b>	<b>No</b>	<b>No</b>	
Supplier to the US Armed Services & VA	<b>Yes</b>	<b>No</b>	<b>No</b>	<b>No</b>	
Number of Supplying Labs	<b>150</b>	12	12	150	Use your favorite lab to build your dreamTAP
PDAC Approved for Medicare	<b>Yes</b>	<b>Yes</b>	<b>No</b>	<b>No</b>	Medicare is about 30-40% of the market. This is a must

Dorsal Fin: Somnodent, Respire, ClearDream, Moses and many more....

Tray Style: EMA, Norval, Pantera, Silent Night